

Pretika®

Hand & Foot Micro-Power Buffer

Soften and Smooth Skin

Instructions
for Use
& Warranty
Information

Model HF111

Pretika®

Merging Science
with Skin Care®

Pretika®

Hand & Foot Micro-Power Buffer

Soften and Smooth Skin

*Congratulations and thank you for being a customer of Pretika!
Before you start, please take a moment to get to know your Hand
and Foot Micro-Smoothing Buffer and how it works. For important
information about safety and how to use the device, read the full
Instructions for Use booklet. Contact the experts at Pretika for any
assistance in helping you get started.*

Customer Service Contact:

1-949-481-8818

Email: customerservice@pretika.com

Or visit our website at www.pretika.com

Hours of Operation:

Monday – Friday, 8:00 AM to 4:00 PM Pacific Standard Time


IMPORTANT INFORMATION

It is our experience that if a product uses batteries, many of the problems thought to be product defects are simply caused by weak batteries. To save any unnecessary inconvenience and expense, we encourage you to try replacement batteries before you return the product.

Micro-Power Buffer Features

To Install or Replace Batteries

1. Make sure the buffer is turned off.
2. Remove back cover by placing thumb on battery cover and slide off with downward push.
3. Carefully pull out batteries (if replacing) and insert 2 new “AA” batteries.
4. Use the markings in the battery compartment as a guide to position the (+) positive or (-) negative
5. End of each battery
6. Replace back cover by sliding up cover into closed position.
7. To avoid leakage, check the batteries regularly for deterioration.
8. Remove exhausted batteries and remove the batteries if the Resurfacing Tool is not to be used for a long time.

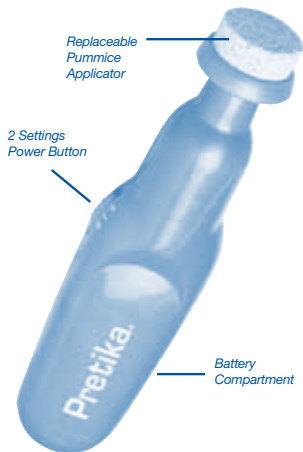

Do not mix old and new batteries. Do not mix alkaline, standard (carbonzinc) or rechargeable (nickel-cadmium) batteries. To dispose of alkaline batteries, place in the trash (normal municipal waste). Recycling of batteries at your local recycling center is good for the environment.

Installing Pumice Stone Head:

Press the Pumice Stone Head (included) onto the connection opening until it locks into position with a click. If replacing the Pumice Stone Head, gently pull Pumice Stone Head to remove.

Daily Treatment

Step 1: Cleanse Skin

Before use, thoroughly cleanse and soak areas to be treated. Rinse well and blot dry.

Step 2: Attach Pumice Stone Applicator. Press applicator to skin and turn on selecting desired speed setting.

Step 3: Use a light circular motion. This step works to remove rough, dry skin, and refine the skin's surface. Work each area for 2-5 minutes, depending on your skin's needs. Wash off residue with warm water and blot dry.

Step 4: For Use in Combination with your Daily Moisturizer Routine

Micro-Power Buffer can be used with many of your daily moisturizing routines. Simply apply your favorite moisturizing after exfoliating and buffing the skin. The Micro-Power Buffer can be used daily, morning and night, but if you have especially sensitive skin or other skin concerns, talk to your doctor or dermatologist about frequency of use.

Caring For Your Resurfacing Tool

1. Wipe handle with slightly damp cloth. Do not immerse Resurfacing Tool in water or any other liquid. Do not open the battery compartment if the appliance is wet to avoid corrosion to the inside parts.
2. Clean attachments after every use. Remove attachments from Resurfacing Tool. Do not allow liquid cleansers or moisturizers to accumulate on attachments or handle.

Replace Pumice Stone Buffer Heads as needed– they can be purchased from your retailer or online at www.pretika.com.

Warning

For external use only. KEEP OUT OF REACH OF CHILDREN. Please be sure to use as directed. As with any product, discontinue use if rash, redness or itching occurs. Consult physician if irritation persists. Use this appliance only for its intended use as described in the Instructions For Use booklet. Do not apply extreme pressure when using unit or when replacing attachments. This could damage product. The attachments are delicate and must be attached and removed with care. Not following operating instructions carefully may result in damage to the unit and attachments. Basic safety precautions should always be followed when using any battery-operated product, especially when children are present. AVOID DIRECT CONTACT WITH EYES AND HAIRLINE/HAIR. IF SKIN IRRITATION OR REDNESS OCCURS, DISCONTINUE USE.

Limited 90 Day Warranty

Your Pretika product is warranted to be free from defect in material and workmanship for a period of 90 days after its original purchase date under normal use. This warranty extends only to the original retailer purchaser with original proof of purchase and only when purchased from an authorized Pretika retailer or reseller. This warranty does not cover products determined by Pretika Corporation to be damaged by any of the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorized persons
- Use with unauthorized accessories or use other than instructed
- Use of incorrect charging adaptor types
- Any other conditions beyond our control

Some states do not allow the exclusion or limitation of incidental or consequential damages or allow limitations on how long an implied warranty lasts, so the above limitations or exclusions may not apply to you.

This warranty gives you specific legal its negligence. Pretika's entire liability for any defective product shall not exceed the purchase price for the defective product. Pretika disclaims any liability for general, consequential, incidental or special damages, even if aware of the likelihood of such damages.

During Warranty Period with Proof of Purchase (store receipt)

Should your product become defective during the warranty period, package the item securely and return it prepaid/insured– along with proof of purchase to the address shown below:

Outside Warranty Period

Should your product become defective outside the warranty period– or you do not have store receipt, package the item securely and return it prepaid/insured along with a check or money order payable to Pretika Corporation in the amount of \$20.00 to cover handling and postage charges. Be sure to include legible name, street address, email and phone number information and send to the address below:

Return to:
Pretika Corporation
Warranty Repair Department
12215 Holly Street
Riverside, CA 92509 USA

To insure prompt repair or replacement, provide complete, legible name, street address, email and phone number information and a note indicating the nature of the product defect. We will Repair or Replace (at our sole discretion) product at no charge and return via ground shipment.

Save these instructions for future reference.



Pretika®
Hand & Foot
Micro-Power Buffer

Soften and Smooth Skin

Pretika Corporation

Riverside, CA 92509

p 949.481.8818 | **f** 949.481.8828

e customerservice@pretika.com

w www.pretika.com

 [Facebook.com/pretikaskincare](https://www.facebook.com/pretikaskincare)

 [Twitter.com/pretikaskincare](https://twitter.com/pretikaskincare)

 www.youtube.com/pretikaskincare